

CASPAR LANT

Intermediate Experimental Physics II

Section: 002

Date Performed: February 2, 2015

Date Due: February 9, 2015

Partner: Sam P. Meier

Professor: Prof. Andrew Kent

Instructor: David Mykytyn

Date: February 8, 2016.

The Objective of this week's experiment was to put our vast theoretical knowledge of lenses to application. It is always a remarkable think to see what was once pure abstraction validated though rigorous scientific experimentation.

1. THEORETICAL BACKGROUND/ ABSTRACT

Just over a century ago, Danish physicist and renowned goalkeeper Niels Bohr postulated the discrete nature of the energy states available to the sole electron in the hydrogen atom.

$$(1) \quad \frac{1}{s_i} + \frac{1}{s_o} = \frac{1}{f}$$

$$(2) \quad \frac{1}{f} = (n - 1) \left(\frac{1}{R_1} - \frac{1}{R_2} \right)$$

2. EXPERIMENTAL PROCEDURE

- (1) Carefully remove the extremely expensive tube from its bubble wrap cocoon.
- (2) Delicately slot the tube into the socket of the supplied cable, and place it in the oven.
- (3) Attach the other end of the cable to the propriety apparatus
- (4) Do not break the tube! It's expensive!
- (5) Turn on the measurement apparatus, as well as the oscilloscope.
- (6) Set the switch to reset.
- (7) Set the value of U_3 to 1.5V.
- (8) Wait about 20 minutes for the Frank-Hertz tube to reach its operating temperature $\vartheta_s = 188^\circ$.
- (9) Turn the largest dial to the sawtooth position, and dial U_1 and U_3 such that the function displayed on the 'scope approximates the curve found on the last page of the instruction packet.
- (10) Set up DataStudio with voltage sensors A and B connected to the oscilloscope inputs.
- (11) Create a graph of voltage B vs. voltage A.
- (12) Send a single sawtooth wave into the DataStudio interface by using the dial on the Frank-Hertz power supply. You should see a more detailed view of the initial oscilloscope trace on your DataStudio display.
- (13) Calculate the distance between two peaks using DataStudio's wonderfully intuitive interface! This should come to about 4.9V; the value postulated by Niels Bohr well before your grandparents were alive.
- (14) You're done! Don't break the tube!

3. GRAPHS AND CALCULATIONS

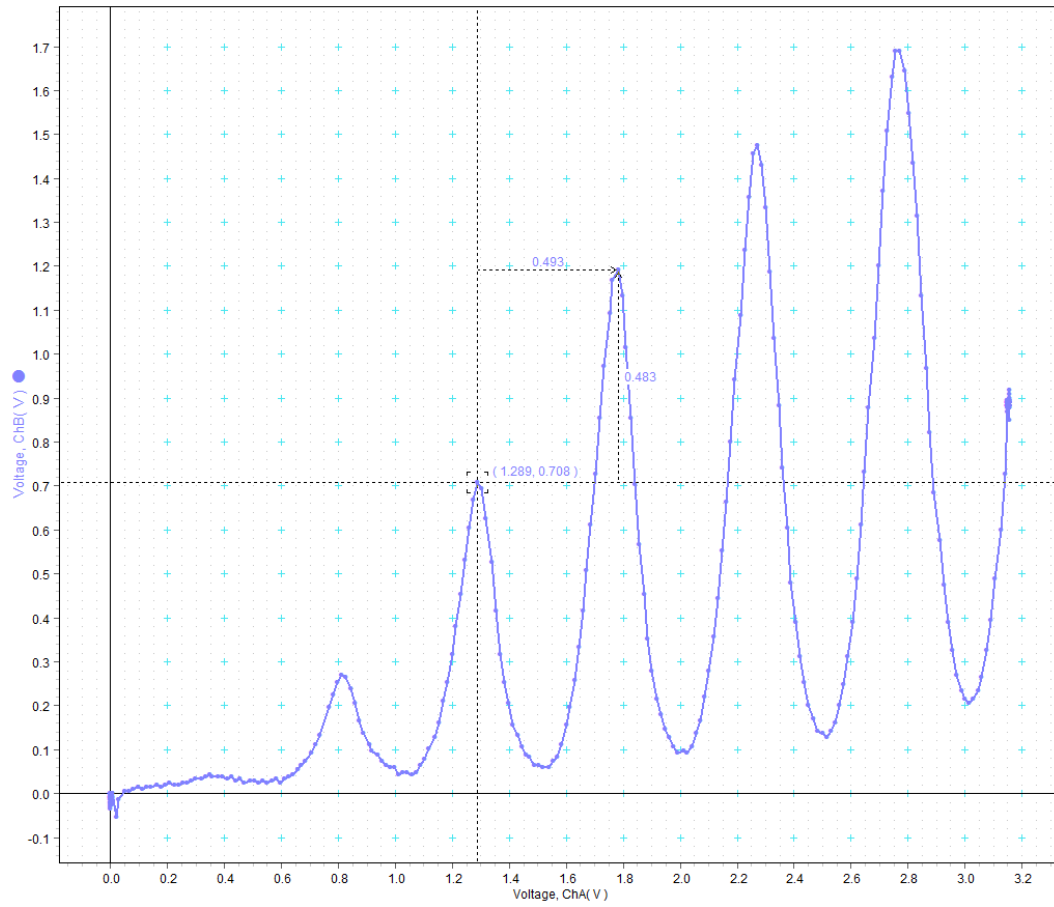


FIGURE 1. Discrete Energy Levels at Intervals of 4.9V

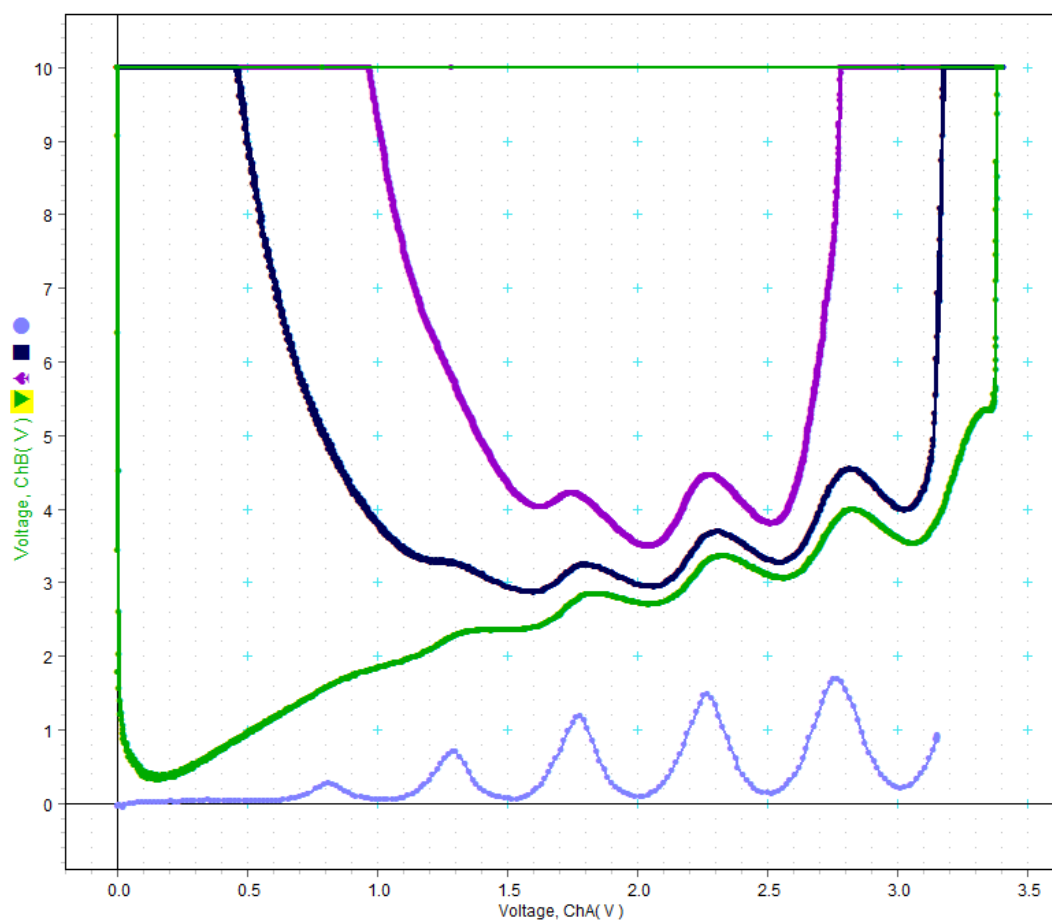


FIGURE 2. Saturation

4. GRAPHS AND CALCULATIONS

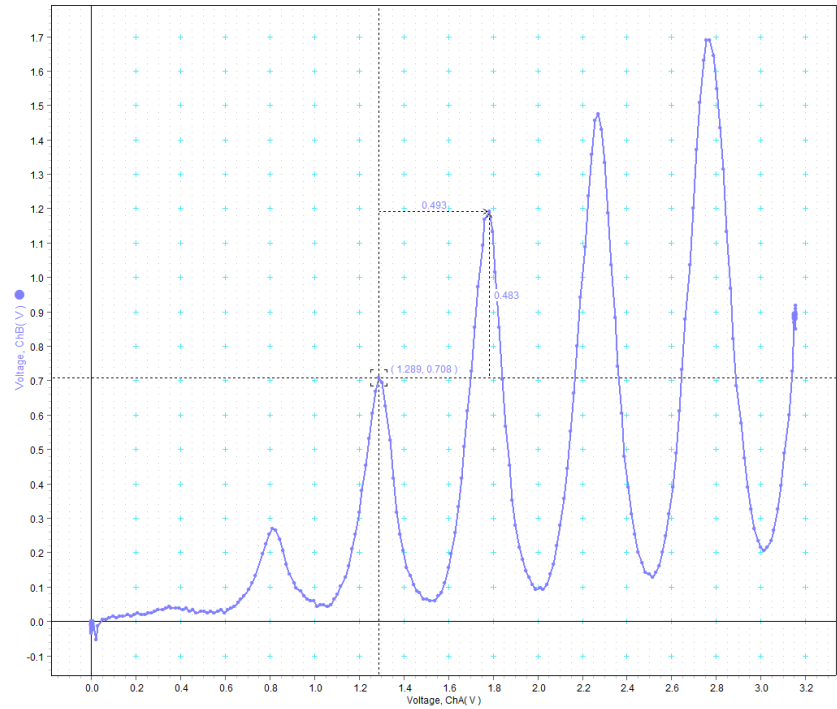


FIGURE 3. Discrete Energy Levels at Intervals of 4.9V

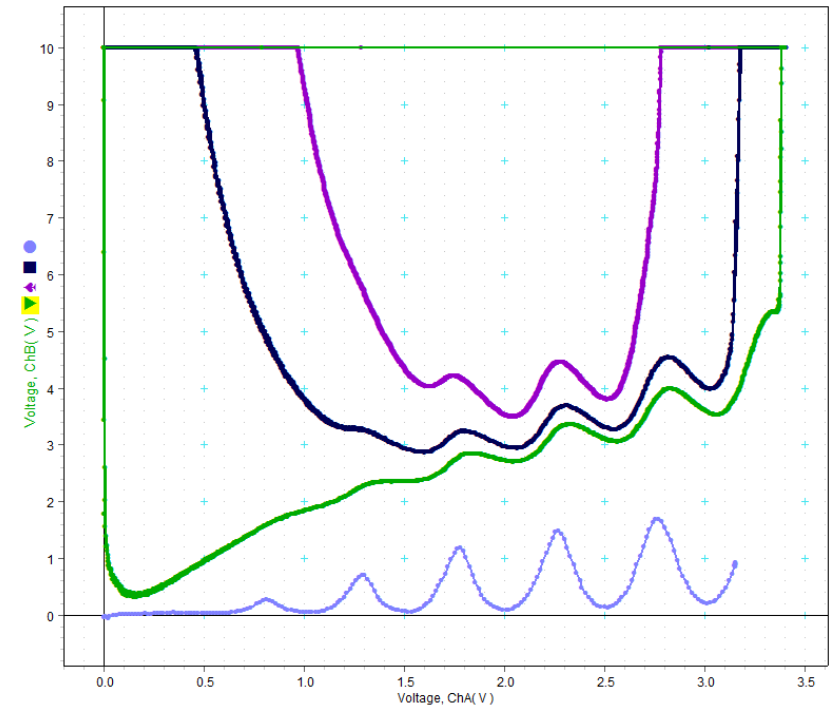


FIGURE 4. Saturation

5. QUESTIONS

- (1) *Why is it better to have the cathode indirectly heated rather than directly heated?*

It is easier to regulate the temperature of the hot cathode if it is heated indirectly. The junction between the heater and the cathode acts as a buffer against fluctuating currents.

- (2) *When the oven temperature is too low, why is there the possibility of a discharge?*

The logarithm of the vapor pressure of a substance is inversely proportional to its temperature. It is important that we run our experiment at sufficiently high temperatures, and low enough corresponding vapor pressure such that vaporized mercury atoms are well-dispersed throughout the tube.

- (3) *Why should the oven temperature not be too high?*

If the oven temperature is too high, we run the risk of damaging the extremely expensive tube at enormous cost. Neglecting glass-breakage, too high a temperature would result in a reduction of the mean-free path of our electrons. This means, physically, that our electrons would experience too many collisions for the experiment to produce any observable effects of quantum energy states.

- (4) *Does the first peak occur at $(U1 + U2) = 4.9$ V? Can you think of reasons as to why it would not?*

I'd say that this engima was due to the presence of the oscilloscope in our circuit. The device drops a constant voltage between the output of the power supply and our lab interface. This has the effect of lowering the voltage at which the first spike occurs, and has no bearing on the distance between consectutive spikes, invariably 4.9V. Revision: I went back into the lab and ran the experiment again, this time with the oscilloscpoe off. The presence of the oscilloscope in my circuit seems only to affect the voltage of the peaks, and not the voltage at which the peaks occur. See the digram below: The data in brown was collected with the oscillospe on, and the data in blue was not.

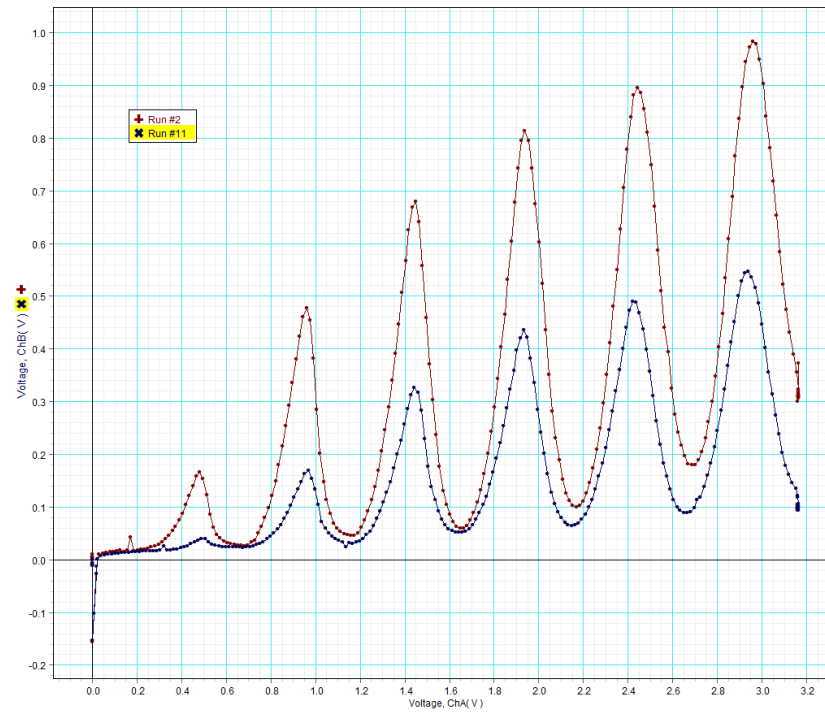


FIGURE 5. NoScope!

You'll This was again confirmed when I probed the oscilloscope terminals with a multimeter with the Frank-Hertz powered down, and found a voltage drop of effectively zero, as shown in the photo below

6. ERROR ANALYSIS